

A First for Drive DeVilbiss, the New...

REMOTE “AUTOFOLD” POWERCHAIR

 **AUTOFOLD**



Following the huge success of our manual folding powerchair, we are delighted to launch our first **AUTOFOLD Powerchair**. This excellent new addition to the range features a remote control function that enables the chair to be easily folded or unfolded in seconds at the simple touch of a button.



LL498

www.drivedevilbiss.co.uk

drive

 **DeVilbiss**
HEALTHCARE

New to Drive DeVilbiss, the AUTOFOLD Powerchair includes several new features to provide a comfortable and easy ride to enhance the user experience. It has a stylish modern design that certainly catches the eye and remote folds to an excellent compact size, ideal for transportation in a car, train, coach or airplane storage.

The AUTOFOLD Powerchair includes Shock Absorbing front suspension to help provide extra comfort for a smoother ride and it can take a very impressive max user weight of 23.6st (150kg).



Codes:

Code	Description
AFPC17BLK	Autofold Powerchair

Specifications:

Unfolded (L x W x H) (cm / in): 101 x 63 x 98 cm / 39.8" x 24.8" x 38.6"	Backrest Height: 48 -58cm (18.8"-22.8")
Folded (L x W x H) (cm / in): 76 x 63 x 42 cm / 30" x 24.8" x 16.5"	Front Wheels: 20 x 5 cm (7.9" x 2")
Seat Height: 51 cm / 20"	Rear Wheels: 30 x 5.7 cm (11.8" x 2.24")
Seat Width x Depth: 44 cm (17.3") x 44 cm (17.3")	Product weight: 28 kg (61.7 lb) excluding battery
Armrest Height: 22cm (8.7")	Battery weight: 1.8 kg (4 lb)

*Speed and range may vary depending upon user weight, type and incline of terrain, weather, battery charge and condition, operating speed and general driving situation.



User friendly
Intelligent controller



Handy underseat
storage



Front wheel
suspension



Flip up armrests



Stylish black/
silver hubs




Remote Fold
Operation

Features:


- » New electric Autofolding function, provides effortless folding and unfolding
- » Remote control fob operation, backed up with a switch on the rear of the powerchair
- » Striking black Aluminium alloy frame, weighing only 28kg (excluding batteries)
- » Front suspension for extra comfort
- » Folds to an excellent compact size, ideal for car, train, coach, and airplane travel
- » Underseat bag as standard for easy to reach storage of your valuables, as well as a backrest pouch for additional items
- » Height adjustable backrest
- » Max Speed 4mph and a maximum range of 7.5miles (12km)*
- » Comfortable and removable seat cushion & padded backrest
- » Anti-Tip wheels and Lap belt for added safety and security
- » Angle adjustable and flip up armrest
- » Flip up footrest (nonslip PU coating)
- » Intelligent & very responsive Joystick controller includes speed and battery gauges, speed control and horn
- » Easy to remove lightweight lithium battery pack
- » On or off board charging
- » Linix brushless DC motors for quality and confidence
- » 24V 10ah LG lithium battery (only requires one battery) with 3-pin charging port
- » Stylish Black Alloy wheels hubs with puncture proof PU tyres
- » Electromagnetic brakes
- » Ground Clearance 6.5cm (2.6")
- » Minimum Turning Radius 100 cm (39")
- » Max climbing angle 6°
- » Max user weight 23.6st (150kg)
- » Spare batteries available providing additional range



Flip up footrest
(nonslip PU coating)

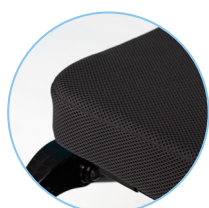

Comfortable padded
seat cushion



If you would like to view this product please contact your local stockist.